




Product Item number	PLUS shower seat 450, manually vertical and horizontal adjustable R7450112035
Color Product description	<div>  </div> <p>Red (RAL 3003)</p> <p>PLUS shower seat 450, with back- and armrests. Manually height adjustable by 7.7" and adjustable horizontally. Armrests are independently height adjustable. The seat and the armrests are foldable. Anthracite cover plate for vertical wall track. Materials: Comfortable polyurethane foam over a plywood core, fitted to an aluminium profile. The non-slip gripping surface of the armrests are made of soft SEBS. For mounting to horizontal wall track (not included).</p> <p>PLUS shower seat 450, with back- and armrests. Manually height adjustable by 7.7" and adjustable horizontally. Armrests are independently height adjustable. The seat and the armrests are foldable. Anthracite cover plate for vertical wall track. Materials: Comfortable polyurethane foam over a plywood core, fitted to an aluminium profile. The non-slip gripping surface of the armrests are made of soft SEBS. For mounting to horizontal wall track (not included).</p>
Adjustability Adjustment height Product capacity	<p>Vertical and horizontal adjustable</p> <p>6,89 inch</p> <p>Maximum user mass (seat): 440 lbs. Maximum load for armrests is 55 kg. When the seat is loaded with 440 lbs, each wall fastener is loaded with 37 lbs.</p> <p><i>The seat has been tested statically with safety factor 1.5, which is a requirement according to the applied assistive products standards.</i></p>
Materials	<p>The shower seat is made of polyurethane (PU) foam with a plywood core and this is covered with a moulded coated surface that is protected using a clear PU lacquer. The backrest is made of polyurethane (PU) foam protected using a clear PU lacquer. The arms are primarily constructed of 3/8" anodized aluminium, with a non-slip gripping surface at the end made of soft SEBS and incorporating a glass-reinforced nylon core. The hinge arm is also made of soft SEBS. Internal components and fasteners are primarily made of stainless steel. The hinge barrel is capped with ABS.</p>

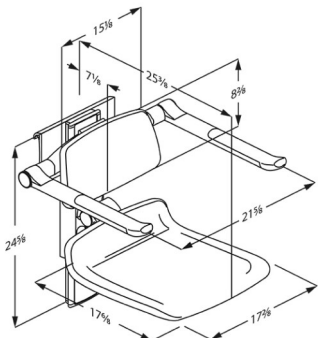

PRESSALIT®

Product Item number	PLUS shower seat 450, manually vertical and horizontal adjustable R7450112035
Ergonomy Features	<p>Ergonomically designed seat - contact surface are 17.7" x 17.2". Soft and slip-resistant material that has a warm feel. The oval shape of the armrests gives an easy and comfortable surface to lean on and push against.</p> <p>The seat, backrest and armrests can be adjusted simultaneously thanks to a safe and easy click system. The armrests can also be adjusted independently of the seat and backrest. The seat and armrests can fold up when not in use or to maximize floor space.</p>
Warranty Cleaning	<p>Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.</p> <p>This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: Erisan OXY+ Tvättsprit 70% CaviCide surface disinfectant decontaminant cleaner Rodalon Desinfektionsspray Terralin protect (SCHÜLKE+) Dismonzon Plus (BODE Chemie GmbH)</p> <p>NON-DISINFECTANTS: Wetrok Caledor (Wetrok) Petrosol Citro (Büfa)</p> <p>Important: Polyurethane surfaces of PLUS shower seats and PLUS backrests contains an antibacterial surface, which will be compromised if cleaned with anything other than soap and water.</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS: Pressalit recommends that water deposits (Calcium, lime etc.) be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in plain water.</p>



PRESSALIT®

Product	PLUS shower seat 450, manually vertical and horizontal adjustable
Item number	R7450112035

pressalit.com

Pressalit A/S
Pressalitvej 1
DK-8680 Ry
Tel.: +45 8788 8788
Fax: +45 8788 8789
E-mail: pressalit@pressalit.com

PRESSALIT®