




# PRESSALIT®

<b>Product</b> <b>Item number</b>	PLUS sink bracket with wired hand control, electrically height manually adjustable and adjustable horizontally R4781-93
<b>Product description</b>	<p>EAN:5708590316185 FDA: ILW/D142473</p>  <p>PLUS sink brackets are available in electrically, pneumatically and manually height adjustable versions. They are suitable for wall-mounting or for fitting to Pressalit PLUS wall tracks.</p> <p>PLUS sink bracket, powered with wired remote control, electrically height adjustable by 11.8" and adjustable horizontally. Recommendation: Flexible feed and waste system R2084/R2085/R2087. Materials: Aluminium and ABS polymer. Incl. mounting bolts for sink. For mounting to horizontal wall track (not included).</p>
<b>Adjustability</b> <b>Adjustment height</b> <b>Product capacity</b>	<p>Vertical and horizontal adjustable</p> <p>11,82 inch</p> <p>Maximum load: 220 lbs. Suitable for sinks weighing no more than 49 lbs. When the sink is loaded with 220 lbs, each wall fastener will be carrying 49 lbs. <i>Sink brackets need to be able to carry the weight of the sink, as well as support the weight of the person who is leaning against the sink or who is using the sink to rise from seated position.</i> <b>Neither sink nor sink bracket are intended to carry the full body weight of a person, and the test methods applicable to sink brackets reflect this: ISO 17966 prescribes a static test on the front edge of the sink at 75% of maximum user weight, i.e. when the maximum load on a sink bracket is given as 220 lbs and the wash basin weighs 49 lbs, it is – according to ISO 17966 – the equivalent of a user weighing 229 lbs leaning on the front edge of the wash basin.</b></p>
<b>Speed up</b> <b>Noise level</b>	<p>Speed up: 0,79 inch/s Speed down: 0,91 inch/s</p> <p>45 dB</p>

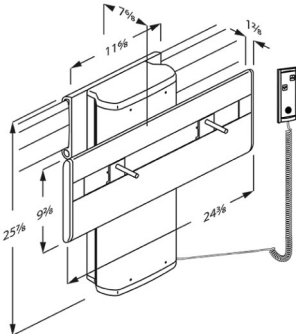

<b>Product</b>	PLUS sink bracket with wired hand control, electrically height manually adjustable and adjustable horizontally
<b>Item number</b>	R4781-93
<b>Materials</b>	The sink bracket is primarily made of 3/8" anodized aluminium. The ends of the bracket cover are capped with ABS polymer. The bracket is available with a lever control consisting of an ABS polymer switch at the end of an anodized aluminium arm
<b>Ergonomy</b>	The height of the sink is easily adjusted the lever control. The sink bracket is compatible with all types of bolt-mounted sink. It's rounded, sleek design makes it easy to clean thoroughly.
<b>Features</b>	The user can easily make height adjustments independently using the electrically operated sink bracket.
<b>Mains connection (V)</b>	120
<b>Actuator (motor)</b>	LA40 (Linak) 4000N, IPX6
<b>Control box</b>	SMPS14002
<b>Wired remote control</b>	HB3 (Linak) IPX6
<b>Power consumption (maximum load)</b>	Power consumption by maximum load: 75 W Power consumption by standby: 1,5 W
<b>Warranty</b>	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.



# PRESSALIT®

<b>Product</b> <b>Item number</b>	PLUS sink bracket with wired hand control, electrically height manually adjustable and adjustable horizontally  R4781-93
<b>Cleaning</b>	<p>This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p><b>DISINFECTANTS:</b> Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70%</p> <p><b>NON-DISINFECTANTS:</b> Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p><b>DECALCIFICATION OF ALUMINUM PARTS:</b> Pressalit recommends that water deposits (Calcium, lime etc.) be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in plain water.</p>



# PRESSALIT®

<b>Product</b>	PLUS sink bracket with wired hand control, electrically height manually adjustable and adjustable horizontally
<b>Item number</b>	R4781-93

---

[pressalit.com](http://pressalit.com)

---

Pressalit A/S  
Pressalitvej 1  
DK-8680 Ry  
Tel.: +45 8788 8788  
Fax: +45 8788 8789  
E-mail: [pressalit@pressalit.com](mailto:pressalit@pressalit.com)

# PRESSALIT®