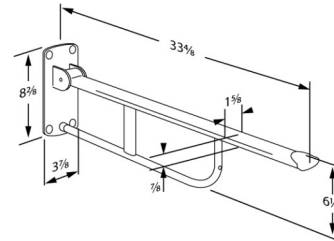




# PRESSALIT®

<b>Product</b> <b>Item number</b>	VALUE fold down grab bar, fixed height R1100000
<b>Color</b>	VVS:615722300 NRF:6191642 RSK:7972410 EAN:5708590176406 HMI: 25411 NBS:N13:436 DB: 2076622 FDA:ILS/D135866 White
<b>Product description</b>	<p>           VALUE fold down grab bar, fixed, foldable, length 33.5".            The ergonomic design gives a secure and comfortable grip as well as a good sense of where the support arm ends.            Incl. fixing bolts for brick wall. Plugs depend on wall type and must be ordered separately.         </p> <p>           VALUE fold down grab bar, fixed, foldable, length 33.5".            The ergonomic design gives a secure and comfortable grip as well as a good sense of where the support arm ends.            Materials: Rilsan-coated steel, pressure-cast aluminium, plastic.            Incl. fixing bolts for brick wall. Plugs depend on wall type and must be ordered separately.         </p>
<b>Adjustability</b>	Fixed mounting
<b>Product capacity</b>	Maximum load: 265 lbs. When the fold down grab bar is loaded with 265 lbs, each mounting bolt is loaded with 534 lbs.
<b>Materials</b>	Basic material: Steel Surface material: Rilsan End cap: Nylon (PA)





<b>Product</b>	VALUE fold down grab bar, fixed height
<b>Item number</b>	R1100000
<b>Ergonomy</b>	<p>Provides firm support when you are sitting, standing, or transferring to and from a toilet.</p> <p>Good, broad support surface for the hand and arm. The rounded form gives a firm and comfortable grip.</p>
<b>Features</b>	<p>The grab bar can fold up when not in use or to make more space.</p>
<b>Warranty</b>	<p>Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.</p>
<b>Cleaning</b>	<p>This product should be cleaned with hot water or a mild cleaning solution. Do not use corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>For surface disinfection, chlorine-based disinfectants may be used according to the product specification and safety precautions for these products.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p><b>DISINFECTANTS:</b> Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% WipeClean UniversalCleaning Chlorine Disinfection Wipes Klorilli (Kiilto Clean)</p> <p><b>NON-DISINFECTANTS:</b> Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p>



# PRESSALIT<sup>®</sup>

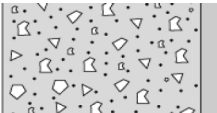



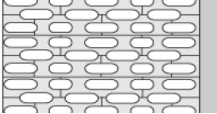
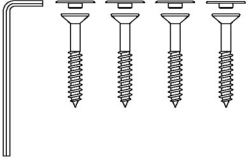
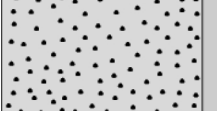
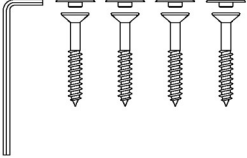
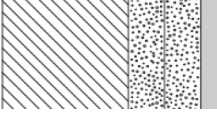
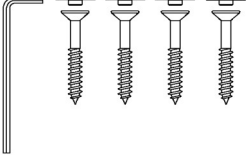
<b>Product</b>	VALUE fold down grab bar, fixed height
<b>Item number</b>	R1100000
<b>Declaration of conformity</b>	This product is in compliance with all relevant requirements of regulation (EU) 2017/745 on medical devices.



PRESSALIT®

<b>Product</b>	VALUE fold down grab bar, fixed height
<b>Item number</b>	R1100000

## Mounting in different wall types

	<p><b>Concrete &gt;B25</b></p> <p>V8593 - 4 countersunk bolts Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.</p>	
	<p><b>Solid Brick &gt;= Mz12</b></p> <p>V8593 - 4 countersunk bolts Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.</p>	
	<p><b>Solid Sand Lime Brick &gt;= Hlz12</b></p> <p>V8593 - 4 countersunk bolts Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.</p>	
	<p><b>Aerated Concrete P4</b></p> <p>V8593 - 4 countersunk bolts Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.</p>	
	<p><b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood</b></p> <p>V8593 - 4 countersunk bolts Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.</p>	



# PRESSALIT®

<b>Product</b> <b>Item number</b>	VALUE fold down grab bar, fixed height R1100000
	<p><b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood</b></p> <p>Wall fasteners (4 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.</p>
	<p><b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (.08 inch)</b></p> <p>Wall fasteners (4 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.</p>
	<p><b>Drywall (2 x 12,5 mm) water-resistant</b></p> <p>Mounting in this wall type is not recommended - risk of fatal injury to persons and walls.</p>

Bolts for wall fixing are included. It is the responsibility of the installer to ensure that suitable plugs are used. Please see mounting instruction.

[pressalit.com](http://pressalit.com)

Pressalit A/S  
 Pressalitvej 1  
 DK-8680 Ry  
 Tel.: +45 8788 8788  
 Fax: +45 8788 8789  
 E-mail: [pressalit@pressalit.com](mailto:pressalit@pressalit.com)

# PRESSALIT®