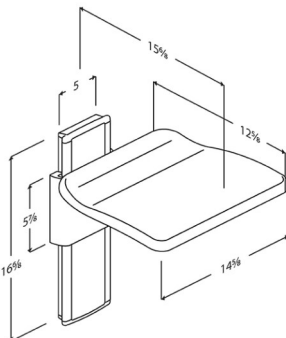




PRESSALIT®

Product Item number	PLUS folding seat, manually height adjustable R5520112000
Color	VVS:777531200 NRF:7501003 RSK:7973927 EAN:5708590306360 HMI: 130158 NBS:N13:448 FDA:INN/D135865 White
Product description	 <p>PLUS folding seat. Manually height adjustable by 9.4". Anthracite cover plate for vertical wall track. The seat is foldable. Excl. mounting kit. The folding seat is not suitable for use in the shower area.</p> <p>PLUS folding seat. Manually height adjustable by 9.4". Anthracite cover plate for vertical wall track. Materials: Aluminium, thermoset plastic, plastic. Excl. mounting kit.</p>
Adjustability Adjustment height Product capacity	Height adjustable 9,46 inch Maximum user mass: 584 lbs. When the seat is loaded with 584 lbs, each wall fastener is loaded with 51 lbs. <i>The seat has been tested statically with safety factor 1.5, which is a requirement according to the applied assistive products standards.</i>
Materials Ergonomy Features	<p>The folding seat is made out of thermoset plastic, and the mounting and suspension bracket are made from 3/8" anodised aluminium. The end caps are made from ABS plastic.</p> <p>Ergonomically shaped seat, which offers good and comfortable support. The seat is suitable for situations where additional seating is required and space is limited, - e.g. in kitchens, lifts and hallways.</p> <p>The seat can be folded away when not in use, - or to free up space. Stays up by itself when folded away.</p>



PRESSALIT®

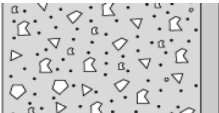
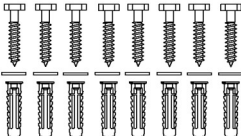

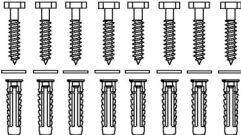
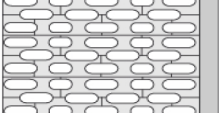
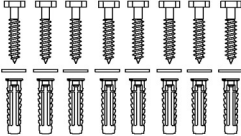
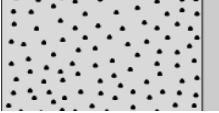
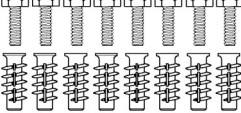
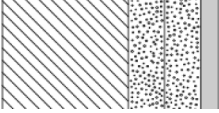
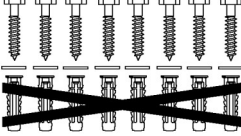
Product	PLUS folding seat, manually height adjustable
Item number	R5520112000
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.
Cleaning	<p>This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: CaviCide</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS: Pressalit recommends that water deposits (Calcium, lime etc.) be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in plain water.</p>
Declaration of conformity	This product is in compliance with all relevant requirements of regulation (EU) 2017/745 on medical devices.



PRESSALIT®

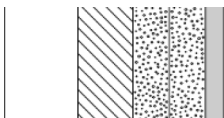
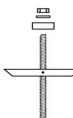


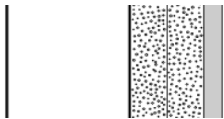
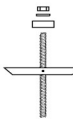


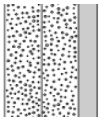
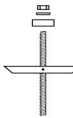


Product	PLUS folding seat, manually height adjustable
Item number	R5520112000

Mounting in different wall types

	<p>Concrete >B25</p> <p>V8610 - Mounting kit, 8 lag bolts Ø8 x 70, plugs SX10. To be ordered separately.</p> 
	<p>Solid Brick >= Mz12</p> <p>V8610 - Mounting kit, 8 lag bolts Ø8 x 70, plugs SX10. To be ordered separately.</p> 
	<p>Solid Sand Lime Brick >= Hlz12</p> <p>V8610 - Mounting kit, 8 lag bolts Ø8 x 70, plugs SX10. To be ordered separately.</p> 
	<p>Aerated Concrete P4</p> <p>V8660 - Mounting kit, 8 countersunk bolts M8 x 35, plugs FTP M 8. To be ordered separately.</p> 
	<p>Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood</p> <p>V8610 - Mounting kit, 8 lag bolts Ø8 x 70, plugs SX10. NB! Plugs are not needed. To be ordered separately.</p> 



PRESSALIT®

Product Item number	PLUS folding seat, manually height adjustable R5520112000		
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood Wall fasteners (8 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.	  	
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (.08 inch) Wall fasteners (8 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.	  	
	Drywall (2 x 12,5 mm) water-resistant Wall fasteners (8 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.	  	

pressalit.com

Pressalit A/S
 Pressalitvej 1
 DK-8680 Ry
 Tel.: +45 8788 8788
 Fax: +45 8788 8789
 E-mail: pressalit@pressalit.com

PRESSALIT®