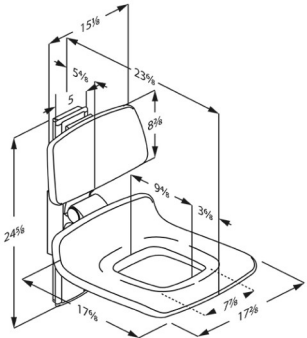




Product Item number	PLUS shower seat 450 with aperture, manually height adjustable R7421112000
Color	VVS:777535600 NRF:6311167 EAN:5708590311692 HMI: 71802 NBS:N13:448 FDA: INN/D135865 White
Product description	 <p>PLUS shower seat 450 with aperture, with backrest. Manually height adjustable by 14.2" . The seat is foldable. Anthracite cover plate for vertical wall track. Materials: Comfortable polyurethane foam over a plywood core, fitted to an aluminium profile. Excl. mounting kit.</p> <p>PLUS shower seat 450 with aperture, with backrest. Manually height adjustable by 14.2" The seat is foldable. Anthracite cover plate for vertical wall track. Materials: Comfortable polyurethane foam over a plywood core, fitted to an aluminium profile. Excl. mounting kit.</p>
Adjustability Adjustment height Product capacity	Height adjustable 14,18 inch Maximum user mass: 440 lbs. When the seat is loaded with 440 lbs, each wall fastener is loaded with 40 lbs. <i>The seat has been tested statically with safety factor 1.5, which is a requirement according to the applied assistive products standards.</i>
Materials Ergonomy	The shower seat is made of polyurethane (PU) foam with a plywood core and this is covered with a moulded coated surface that is protected using a clear PU lacquer. The backrest is made of polyurethane (PU) foam protected using a clear PU lacquer. Internal components and fasteners are primarily made of stainless steel. The hinge barrel is capped with ABS. Ergonomically designed seat - contact surface are 17.7" x 17.2". Soft and slip-resistant material that has a warm feel.



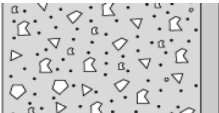
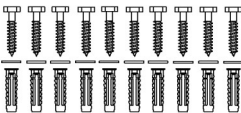

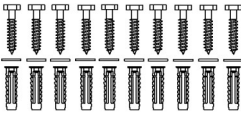
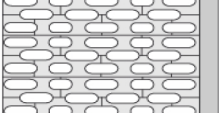
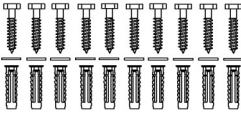
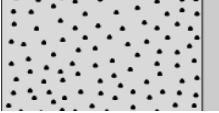
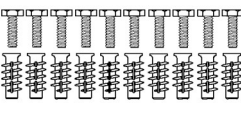
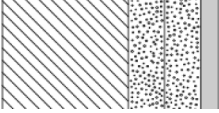
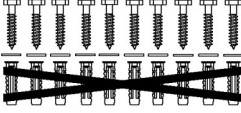
Product	PLUS shower seat 450 with aperture, manually height adjustable
Item number	R7421112000
Features	The seat can fold up when not in use or to maximize floor space. The seat and backrest can be adjusted simultaneously thanks to a safe and easy click system.
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.
Cleaning	<p>This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: Erisan OXY+ Tvättsprit 70% CaviCide surface disinfectant decontaminant cleaner Rodalon Desinfektionsspray Terralin protect (SCHÜLKE+) Dismonzon Plus (BODE Chemie GmbH)</p> <p>NON-DISINFECTANTS: Wetrok Caledor (Wetrok) Petrosol Citro (Büfa)</p> <p>Important: Polyurethane surfaces of PLUS shower seats and PLUS backrests contains an antibacterial surface, which will be compromised if cleaned with anything other than soap and water.</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS: Pressalit recommends that water deposits (Calcium, lime etc.) be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in plain water.</p>



PRESSALIT®

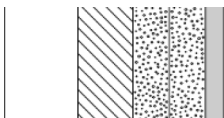
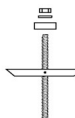

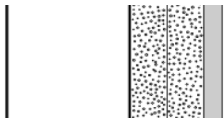
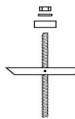

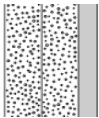
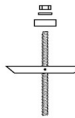

Product	PLUS shower seat 450 with aperture, manually height adjustable
Item number	R7421112000

Mounting in different wall types

	<p>Concrete >B25</p> <p>V8612 - Mounting kit, 10 lag bolts Ø8 x 70, plugs SX10. To be ordered separately.</p> 
	<p>Solid Brick >= Mz12</p> <p>V8612 - Mounting kit, 10 lag bolts Ø8 x 70, plugs SX10. To be ordered separately.</p> 
	<p>Solid Sand Lime Brick >= Hlz12</p> <p>V8612 - Mounting kit, 10 lag bolts Ø8 x 70, plugs SX10. To be ordered separately.</p> 
	<p>Aerated Concrete P4</p> <p>V8662 - Mounting kit, 10 countersunk bolts M8 x 35, plugs FTP M 8. To be ordered separately.</p> 
	<p>Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood</p> <p>V8612 - Mounting kit, 10 lag bolts Ø8 x 70, plugs SX10. NB! Plugs are not needed. To be ordered separately.</p> 



PRESSALIT®

Product Item number	PLUS shower seat 450 with aperture, manually height adjustable R7421112000		
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood Wall fasteners (10 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.		
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (.08 inch) Wall fasteners (10 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.		
	Drywall (2 x 12,5 mm) water-resistant Wall fasteners (10 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.		

pressalit.com

Pressalit A/S
 Pressalitvej 1
 DK-8680 Ry
 Tel.: +45 8788 8788
 Fax: +45 8788 8789
 E-mail: pressalit@pressalit.com

PRESSALIT®