




<b>Product</b> <b>Item number</b>	PLUS shower seat 450 with aperture, manually vertical and horizontal adjustable R7441182000
<b>Color</b>          <b>Product description</b>	<div> <div>White</div>  </div> <p>PLUS shower seat 450 with aperture, with backrest. Manually height adjustable by 14.7" and adjustable horizontally. The seat is foldable. Pale gray cover plate for vertical wall track.          Materials: Comfortable polyurethane foam over a plywood core, fitted to an aluminium profile.          For mounting to horizontal wall track (not included).</p> <p>PLUS shower seat 450 with aperture, with backrest. Manually height adjustable by 14.2" and adjustable horizontally. The seat is foldable. Pale gray cover plate for vertical wall track.          Materials: Comfortable polyurethane foam over a plywood core, fitted to an aluminium profile.          For mounting to horizontal wall track (not included).</p>
<b>Adjustability</b>  <b>Adjustment height</b>   <b>Product capacity</b>	Vertical and horizontal adjustable 14,18 inch   Maximum user mass: 440 lbs. When the seat is loaded with 440 lbs, each wall fastener is loaded with 37 lbs.  <i>The seat has been tested statically with safety factor 1.5, which is a requirement according to the applied assistive products standards.</i>
<b>Materials</b>          <b>Ergonomy</b>	<p>The shower seat is made of polyurethane (PU) foam with a plywood core and this is covered with a moulded coated surface that is protected using a clear PU lacquer. The backrest is made of polyurethane (PU) foam protected using a clear PU lacquer. Internal components and fasteners are primarily made of stainless steel. The hinge barrel is capped with ABS.</p> <p>Ergonomically designed seat - contact surface are 17.7" x 17.2". Soft and slip-resistant material that has a warm feel.</p>

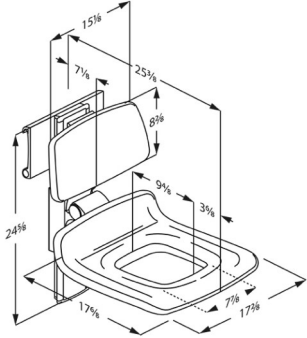

# PRESSALIT®

<b>Product</b> <b>Item number</b>	PLUS shower seat 450 with aperture, manually vertical and horizontal adjustable R7441182000
<b>Features</b>	The seat can fold up when not in use or to maximize floor space. The seat and backrest can be adjusted simultaneously thanks to a safe and easy click system.
<b>Warranty</b>  <b>Cleaning</b>	<p>Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.</p> <p>This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p><b>DISINFECTANTS:</b> CaviCide</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p><b>DECALCIFICATION OF ALUMINUM PARTS:</b> Pressalit recommends that water deposits (Calcium, lime etc.) be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in plain water.</p>

[pressalit.com](http://pressalit.com)

## PRESSALIT®