




PRESSALIT®

Product Item number	PLUS backrest, fixed height R6510000
Color	VVS:777538200 NRF:6190615 EAN:5708590314303 HMI: 71810 NBS:N13:426 FDA:INN/D135865 White
Product description	 <p>PLUS backrest 15.2". For fixed mounting. Incl. mounting kit for brick wall.</p> <p>PLUS backrest 15.2". For fixed mounting. Incl. mounting kit for brick wall.</p>
Adjustability	Fixed mounting
Product capacity	Maximum load: 220 lbs.
Materials	The backrest is made of polyurethane (PU) foam which is covered with a moulded coated surface that is protected using a clear PU lacquer.
Ergonomy	Soft and slip-resistant material that has a warm feel.
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.

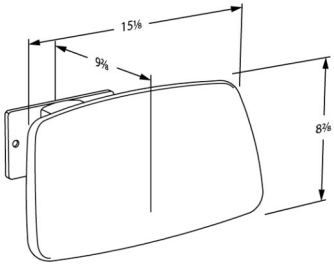

PRESSALIT®

Product	PLUS backrest, fixed height
Item number	R6510000
Cleaning	<p>This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: CaviCide</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS: Pressalit recommends that water deposits (Calcium, lime etc.) be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in plain water.</p>

pressalit.com

PRESSALIT®