




PRESSALIT®

Product Item number	PLUS backrest, vertical and horizontal adjustable R6530112000
Color	VVS:777538600 NRF:6190622 EAN:5708590314488 HMI: 71808 NBS:N13:426 FDA:INN/D135865 White
Product description	 <p>PLUS backrest 15.2". Manually height adjustable by 5.9" and adjustable horizontally. Anthracite cover plate for vertical wall track. For mounting to horizontal wall track (not included).</p> <p>PLUS backrest 15.2". Manually height adjustable by 5.9" and adjustable horizontally. Anthracite cover plate for vertical wall track. For mounting to horizontal wall track (not included).</p>
Adjustability Adjustment height Product capacity	Vertical and horizontal adjustable 5,91 inch Maximum load: 220 lbs.
Materials Ergonomy	The backrest is made of polyurethane (PU) foam which is covered with a moulded coated surface that is protected using a clear PU lacquer. Soft and slip-resistant material that has a warm feel.
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.

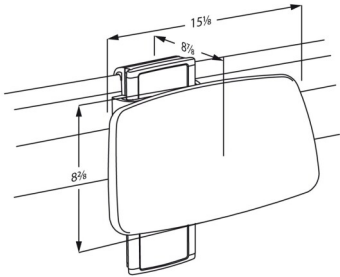

Product	PLUS backrest, vertical and horizontal adjustable
Item number	R6530112000
Cleaning	<p>This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: Erisan OXY+ Tvättspit 70% CaviCide surface disinfectant decontaminant cleaner Rodalon Desinfektionsspray Terralin protect (SCHÜLKE+) Dismonzon Plus (BODE Chemie GmbH)</p> <p>NON-DISINFECTANTS: Wetrok Caledor (Wetrok) Petrosol Citro (Büfa)</p> <p>Important: Polyurethane surfaces of PLUS shower seats and PLUS backrests contains an antibacterial surface, which will be compromised if cleaned with anything other than soap and water.</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS: Pressalit recommends that water deposits (Calcium, lime etc.) be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in plain water.</p>