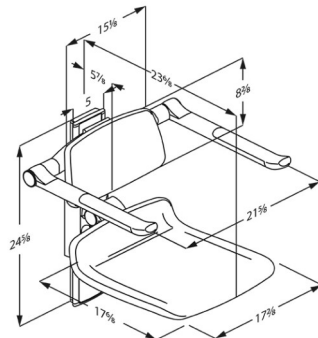




<b>Product</b>	PLUS shower seat 450, manually height adjustable
<b>Item number</b>	R7430112000
<b>Color</b>	<p>VVS: 777535800 NRF: 6311172 EAN: 5708590311876 HMI: 71800 NBS: N13:448 FDA: INN/D135865 White</p> 
<b>Product description</b>	<p>PLUS shower seat 450, with back- and armrests. Manually height adjustable by 6.9". Armrests are independently height adjustable. The seat and the armrests are foldable. Anthracite cover plate for vertical wall track.</p> <p>Materials: Comfortable polyurethane foam over a plywood core, fitted to an aluminium profile. The non-slip gripping surface of the armrests are made of soft SEBS incorporating a glass-reinforced nylon core.</p> <p>Excl. fixing kit.</p> <p>PLUS shower seat 450, with back- and armrests. Manually height adjustable by 6.9". Armrests are independently height adjustable. The seat and the armrests are foldable. Anthracite cover plate for vertical wall track.</p> <p>Materials: Comfortable polyurethane foam over a plywood core, fitted to an aluminium profile. The non-slip gripping surface of the armrests are made of soft SEBS incorporating a glass-reinforced nylon core.</p> <p>Excl. fixing kit.</p>
<b>Adjustability</b>	Height adjustable
<b>Adjustment height</b>	6,89 inch
<b>Product capacity</b>	<p>Maximum user mass (seat): 440 lbs. Maximum load for armrests is 55 kg. When the seat is loaded with 440 lbs, each wall fastener is loaded with 40 lbs.</p> <p><i>The seat has been tested statically with safety factor 1.5, which is a requirement according to the applied assistive products standards.</i></p>



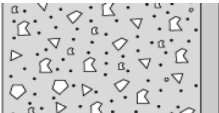
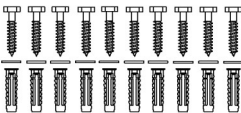

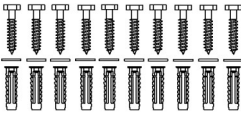
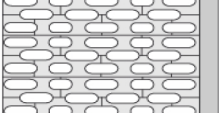
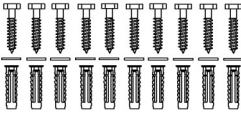
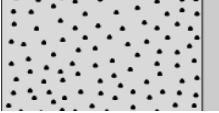
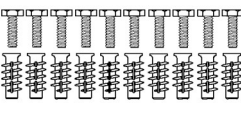
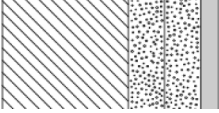
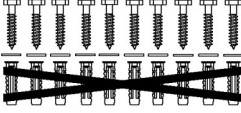
<b>Product</b> <b>Item number</b>	PLUS shower seat 450, manually height adjustable R7430112000
<b>Materials</b>  <b>Ergonomy</b>  <b>Features</b>	<p>The shower seat is made of polyurethane (PU) foam with a plywood core and this is covered with a moulded coated surface that is protected using a clear PU lacquer. The backrest is made of polyurethane (PU) foam protected using a clear PU lacquer. The arms are primarily constructed of 3/8" anodized aluminium, with a non-slip gripping surface at the end made of soft SEBS and incorporating a glass-reinforced nylon core. The hinge arm is also made of soft SEBS. Internal components and fasteners are primarily made of stainless steel. The hinge barrel is capped with ABS.</p> <p>Ergonomically designed seat - contact surface are 17.7" x 17.2". Soft and slip-resistant material that has a warm feel. The oval shape of the armrests gives an easy and comfortable surface to lean on and push against.</p> <p>The seat, backrest and armrests can be adjusted simultaneously thanks to a safe and easy click system. The armrests can also be adjusted independently of the seat and backrest. The seat and armrests can fold up when not in use or to maximize floor space.</p>
<b>Warranty</b>  <b>Cleaning</b>	<p>Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.</p> <p>This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: CaviCide</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS: Pressalit recommends that water deposits (Calcium, lime etc.) be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in plain water.</p>



PRESSALIT®

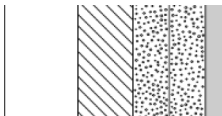
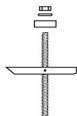

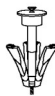
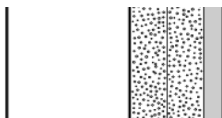
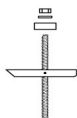

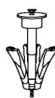
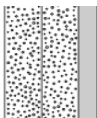
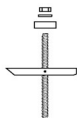

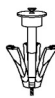
<b>Product</b>	PLUS shower seat 450, manually height adjustable
<b>Item number</b>	R7430112000

## Mounting in different wall types

	<b>Concrete &gt;B25</b> V8612 - Mounting kit, 10 lag bolts Ø8 x 70, plugs SX10. To be ordered separately.	
	<b>Solid Brick &gt;= Mz12</b> V8612 - Mounting kit, 10 lag bolts Ø8 x 70, plugs SX10. To be ordered separately.	
	<b>Solid Sand Lime Brick &gt;= Hlz12</b> V8612 - Mounting kit, 10 lag bolts Ø8 x 70, plugs SX10. To be ordered separately.	
	<b>Aerated Concrete P4</b> V8662 - Mounting kit, 10 countersunk bolts M8 x 35, plugs FTP M 8. To be ordered separately.	
	<b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood</b> V8612 - Mounting kit, 10 lag bolts Ø8 x 70, plugs SX10. NB! Plugs are not needed. To be ordered separately.	



# PRESSALIT®

<b>Product</b> <b>Item number</b>	PLUS shower seat 450, manually height adjustable R7430112000		
	<b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood</b> Wall fasteners (10 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.	  	
	<b>Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (.08 inch)</b> Wall fasteners (10 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.	  	
	<b>Drywall (2 x 12,5 mm) water-resistant</b> Wall fasteners (10 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.	  	

pressalit.com

Pressalit A/S  
 Pressalitvej 1  
 DK-8680 Ry  
 Tel.: +45 8788 8788  
 Fax: +45 8788 8789  
 E-mail: [pressalit@pressalit.com](mailto:pressalit@pressalit.com)

# PRESSALIT®