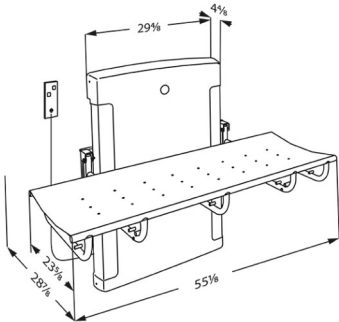




PRESSALIT®

| | |
|--|---|
| Product Item number | Adult changing table 1000, electrically height adjustable R8400000-93 |
| Color | EAN:5708590317939 FDA:INQ/D124319 White |
| Product description |  <p>Discontinued per 31 Dec 2018 and replaced by R8402 Adult changing table 1000, powered, height adjustable 15,8"; from 18.9" - 34.7", with wired hand control Foldable. Outer measurements 28.9" x 55.1". Excl. mounting kit. See Technical specifications (pdf) for information on mounting kits suitable for various wall types.</p> <p>Discontinued per 31 Dec 2018 and replaced by R8402 Adult changing table 1000, powered, height adjustable 15,8"; from 18.9" - 34.7", with wired hand control Foldable. Outer measurements 28.9" x 55.1". Excl. mounting kit. See Technical specifications (pdf) for information on mounting kits suitable for various wall types.</p> |
| Adjustability Adjustment height Product capacity Speed up Noise level | Height adjustable 15,76 inch Maximum load: 330 lbs. When the table is loaded with 330 lbs, each mounting bolt supports 44 lbs. <i>This product has been tested statically with a load of 830 lbs as the IEC 60601-standard requires testing with a safety factor.</i> Speed up: 0,59 inch/s Speed down: 0,79 inch/s 48 dB |



| | |
|---|--|
| Product | Adult changing table 1000, electrically height adjustable |
| Item number | R8400000-93 |
| Materials | The table is made from practical and durable PVC canvas which is phthalate free. The frame is powder-coated steel. Material: Steel, PVC. |
| Ergonomy | Pressalit offers a range of high quality, fixed or power adjusted changing tables designed to provide a stable platform suitable for changing adults and children in a variety of settings: schools, day care centers, institutions, hospitals and private homes. Caregivers can easily set the height of the table, enabling them to obtain the correct working position. This ensures safe showering and changing for the user as well as minimizing the risk of back injury for the caregiver. |
| Features | Electric with wired hand control – the red buttons clearly indicate “Up” and “Down”. The wired hand control comes with a holder for fitting to the wall. The integrated pneumatic cylinders enable the table to be folded out of the way when not in use to maximize floor space in front of the unit. Splash curtain (optional) for the changing table protects the caregiver against splashing. |
| Mains connection (V) | 120 |
| Actuator (motor) | LA31 Careline, (Linak), 4000 N, IPX6 |
| Control box | CB9 Careline, (Linak), IPX6 |
| Wired remote control | HB3 (Linak), IP66 |
| Power consumption (maximum load) | Power consumption by maximum load: 70 W Power consumption by standby: 2 W |
| Warranty | Pressalit makes a commitment to remedy material and manufacturing defects for a period of three years from the invoicing date. NOTICE If the product is placed in an atmosphere heavy in salt or chloride, the warranty will be reduced from the normal 3 years to 1 year. Salt and chloride have a corrosive effect on the products components and will cause a decrease in the life expectancy of the product. |



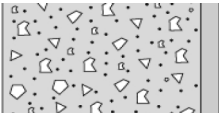
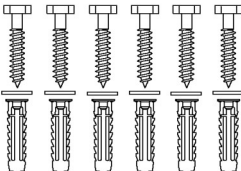

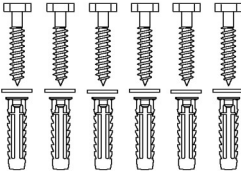
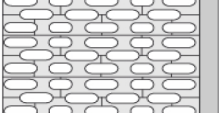
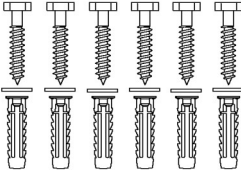
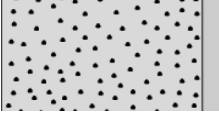
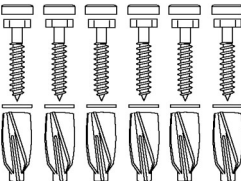
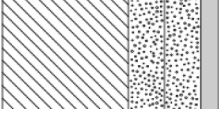
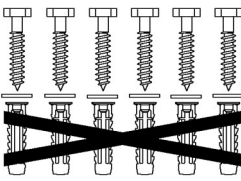
| | |
|--------------------|---|
| Product | Adult changing table 1000, electrically height adjustable |
| Item number | R8400000-93 |
| Cleaning | <p>This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% Accel PREvention Wipes Lingettes</p> <p>NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> |



PRESSALIT®

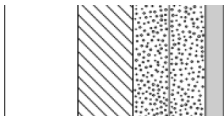
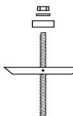

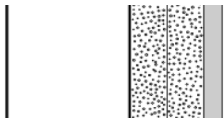
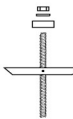

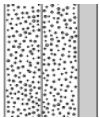

| | |
|--------------------|---|
| Product | Adult changing table 1000, electrically height adjustable |
| Item number | R8400000-93 |

Mounting in different wall types

| | | |
|---|--|---|
|  | <p>Concrete >B25</p> <p>V8608 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs SX12. To be ordered separately.</p> |  |
|  | <p>Solid Brick >= Mz12</p> <p>V8608 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs SX12. To be ordered separately.</p> |  |
|  | <p>Solid Sand Lime Brick >= Hlz12</p> <p>V8608 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs SX12. To be ordered separately.</p> |  |
|  | <p>Aerated Concrete P4</p> <p>V8658 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs GB14. To be ordered separately.</p> |  |
|  | <p>Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood</p> <p>V8608 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs SX12. NB! Plugs are not needed. To be ordered separately.</p> |  |



PRESSALIT®

| | | | |
|---|--|---|---|
| Product Item number | Adult changing table 1000, electrically height adjustable R8400000-93 | | |
|  | Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood Wall fasteners (6 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar. |  |  |
|  | Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (.08 inch) Wall fasteners (6 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar. |  |  |
|  | Drywall (2 x 12,5 mm) water-resistant Mounting in this wall type is not recommended - risk of fatal injury to persons and walls. |  | |

It is the responsibility of the installer to ensure that suitable wall fasteners are used.

pressalit.com