




PRESSALIT®

Product Item number	PLUS replacement shower seat 450 with aperture, manually height adjustable R7401000
Color	White
Product description	<p>VVS: 777534200 EAN: 5708590328430</p>  <p>PLUS shower seat 450 with aperture. Manually height adjustable by 9.1". The seat is foldable. Materials: Comfortable polyurethane foam over a plywood core, fitted to an aluminium profile. For mounting on vertical wall track (not included).</p> <p>PLUS shower seat 450 with aperture. Manually height adjustable by 9.1". The seat is foldable. Materials: Comfortable polyurethane foam over a plywood core, fitted to an aluminium profile. For mounting on vertical wall track (not included).</p>
Adjustability Adjustment height Product capacity	<p>Height adjustable</p> <p>9,06 inch</p> <p>Maximum user mass: 440 lbs. When the seat is loaded with 440 lbs, each wall fastener is loaded with 77 lbs.</p> <p><i>The seat has been tested statically with safety factor 1.5, which is a requirement according to the applied assistive products standards.</i></p>
Materials Ergonomy Features	<p>The shower seat is made of polyurethane (PU) foam with a plywood core and this is covered with a moulded coated surface that is protected using a clear PU lacquer. Internal components and fasteners are primarily made of stainless steel. The hinge barrel is capped with ABS.</p> <p>Ergonomically designed seat - contact surface are 12.2" x 14.4". Soft and slip-resistant material that has a warm feel.</p> <p>The seat can fold up when not in use or to maximise floor space.</p>

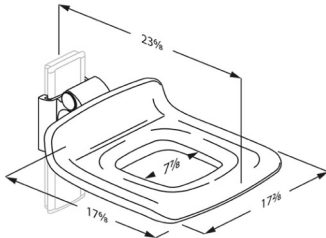

PRESSALIT®

Product	PLUS replacement shower seat 450 with aperture, manually height adjustable
Item number	R7401000
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.
Cleaning	<p>This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: CaviCide</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS: Pressalit recommends that water deposits (Calcium, lime etc.) be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in plain water.</p>

pressalit.com

PRESSALIT®